


Aidacare FL250 Floorline Bed

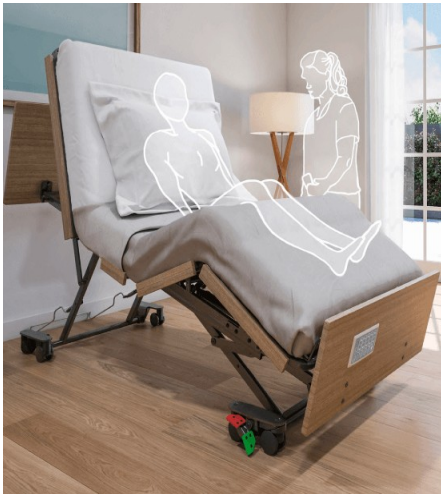

- Under Bed Clearance - While the FL250 achieves a low height of 4', full clearance under the bed is still maintained, minimising the risk of foot or cord impingement. .
- SmartDescent™ Automatically slows descent speed over the final 300mm of lowering to the floor level, allowing possible impingements to be identified and moved.
- Hidden Bed Mechanisms: The bed mechanisms are attached under the platform and are not visible. This helps the room feel less clinical.
- Low Height 4.3': Adding a 5.9' mattress gives you a total height below 11.8' thus significantly reducing out-of-bed fall injury risks.
- Auto-Regression on back & knee rests: Reduces skin shear forces during operation thus improving clinical outcomes and reducing the need for expensive mattress solutions.
- High SWL 39 st & base-extension to 4' & 8': The same bed platform can meet a wider profile of residents.
- BrakeFirst™ Automatic Safety Brake Activation at lowest height: Reduces the risk of injury during bed exit in the morning as the bed will lock when at its lowest height.
- Single Button Seat Position: Easy transfer of weight from lower back to bottom for pressure-care and seated positioning also helps improves lung function for resident.
- ExitSafe™ Single Button Ingress/Exit height: Your staff don't have to think about what height to set the bed. Via one button the vast majority of needs are met.
- Attendant Control Panel - Allows the caregiver to lock out certain functions and also ensures bed operation is not just restricted to the handset.
- Grey Colour Scheme & Narrow-profile frame: Keeps the focus on bedding and furnishings/finishings rather than the bed-frame.
- Underbed light standard: Helps to reduce incidence of resident and caregiver

falls.

- Rubber & Nylon castors: Castors remain on the ground during braking, thus not damaging floors.
- King Single Extension Kit (optional retrofit): Allows you to change the profile of your beds (with minimal cost) to meet resident profile changes.
- KTL coated anti-corrosion steel: Improved durability for lower lifecycle cost.
- Low-maintenance Braking System: Fewer headaches for your MO's!

Warranty Information

- 10 Years Frame
- 3 Years Electrics
- 1 Year Handset

Code	Display Name	Total Length	Total Width	Safe Working Load
UK-BEB046100	Bed - Aidacare FL250 Floorline - Single - 4 Function - 39 st (Bed ends sep) - Incl: UBL & Extension	6'6"	35"	39.5 st
UK-BEB046150	Bed - Aidacare FL250 Floorline - King Single - 4 Function - 39 st (Bed ends sep) - Incl: UBL & Extension	6'6"	3'6"	39.5 st